

THEATER MEDICAL INFORMATION PROGRAM-JOINT (TMIP-J) DATA MIGRATION GUIDE VERSION 1.1

Submitted in Response to:

D/SIDDOMS III Contract No. W74V8H-04-D-0036
Delivery Order 0016
Theater Medical Information Program - Joint (TMIP-J)
Integration Development
Deliverable Item 28

For:

Defense Health Information Management System (DHIMS)
Lt Col Michael Holmes, USAF
5109 Leesburg Pike, Suite 100
Falls Church, VA 22041
(703) 998-6900, ext. 1103

By:

Science Applications International Corporation
Health Solutions Business Unit
5203 Leesburg Pike, Suite 1300, Falls Church, VA 22041
Beth Le Duc, TMIP-J Program Manager, Phone: (703) 824-5808
Amy Ward, TMIP-J Deputy Program Manager, Phone: (703) 664-4834

REVISION HISTORY

Date	Document #	Version	Description	Author
07/24/08	D3-TMIP-5019	0.1	Draft B2 FY07 SP1 delivery to PMO	SAIC Training
08/15/08	D3-TMIP-5019	1.0	Updates for delivery	SAIC Training
09/20/08	D3-TMIP-5019A	1.1	Updates for delivery	SAIC Training

TABLE OF CONTENTS

Sec	tion	l	Page
1.0		SCOPE	1
	1.1	Identification	1
	1.2	Data Migration Overvi	ew1
	1.3	Document Overview	1
2.0		PRE-MIGRATION SE	TUP2
	2.1	Data To Be Migrated	2
	2.2	B1 Data Transfer Via t	he Stem2
		2.2.1 B1 STEM Proce	dure
		2.2.2 Export STEM V	olume
		-	olume4
	2.3	1 0	5
			ystem
			Oracle Version 6 ration Setup 6
	2.4		6
	∠.∓		chine
		-	ration Setup
	2.5		ware Support7
3.0		= -	EXPORT8
4.0		DATA MIGRATION –	IMPORT27
			S/ABBREVIATIONS A-1
AP	PEN	NDIX B: PRE-MIGRA	ATION CHECKLISTB-1
			LIST OF FIGURES
_	ure ure	2-1 STEM Application	Page Screen 3
Fig	ure	3-1 Data Migration set	ıp.exe File Screen9
Fig	ure	3-2 .NET Framework 1	.1 Setup License Agreement Screen 10
Fig	ure	3-3 Installation of Micr	cosoft .NET Framework 1.1 Is Complete Screen11
Fig	ure	3-4 Welcome to the Dat	ta Migration Wizard Screen11
Fig	ure	3-5 Data Migration Gu	ide Dialog Box12
Fig	ure	3-6 Data Migration Re	mote Users Warning Screen 12
Fig	ure	3-7 Choose Working A	nd TMIP B2 Data Container Folders Screen13

Figure 3-8 Choose Working And TMIP B2 Data Container Folders - Completed Screen 14
Figure 3-9 Data Migration Drive Selection Screen
Figure 3-10 Data Migration User Input Screen
Figure 3-11 Select the Product to Migrate Screen
Figure 3-12 List of Components to Migrate Screen
Figure 3-13 Data Migration Setup Status Screen
Figure 3-14 Data Migration Estimated Time to Complete Screen
Figure 3-15 CHCSII-Theater Data Migration Manager Screen
Figure 3-16 Oracle DOS Export Data Screen
Figure 3-17 CHCS II-Theater Process Completed Successfully Screen
Figure 3-18 Data Migration Setup Status Screen
Figure 3-19 Data Migration Setup Status Screen - Continued
Figure 3-20 Task Status Screen
Figure 3-21 Migration Completed Successfully Screen
Figure 3-22 Results of Migration Screen
Figure 3-23 Migration Complete Screen
Figure 3-24 Readme.txt Notepad Screen
Figure 4-1 Data Migration setup.exe Screen
Figure 4-2 Welcome to the Data Migration Wizard Screen
Figure 4-3 Data Migration Error Screen
Figure 4-4 Data Migration Guide Dialog Screen
Figure 4-5 Data Migration Remote Users Warning Screen
Figure 4-6 Choose Working Folder, MyDocs Restore Folder And TMIP B2 Screen 32
Figure 4-7 Select TMIPB2DataContainer Folder from Select File – Completed Screen 33
Figure 4-8 Select the Product to Migrate Screen
Figure 4-9 Select the Product to Migrate Screen
Figure 4-10 List of Components to Migrate Screen
Figure 4-11 Data Migration User Input Screen
Figure 4-12 Task Status Screen
Figure 4-13 Task Status Screen - Continued
Figure 4-14 Setup Status Screen
Figure 4-15 There Is Existing Data in the Application Database Warning Dialog Box 40
Figure 4-16 Data Migration Estimated Time to Complete Screen

Figure 4-17 AHLTA-Theater Data Migration Manager Screen	41
Figure 4-18 DOS Screen for Oracle Product Importing	42
Figure 4-19 AHLTA-Theater Data Migration Manager Screen	43
Figure 4-20 AHLTA-Theater Process Completed Successfully Screen	44
Figure 4-21 Setup Status Continuation Process Screen	45
Figure 4-22 Setup Status My Documents Continuation Process Screen - Continued	45
Figure 4-23 Completed Migration of My Documents Screen - Continued	46
Figure 4-24 DOS Compile Screen	46
Figure 4-25 BoCSharpCompiler Logon Screen	47
Figure 4-26 BoCSharpCompiler Package Screen	48
Figure 4-27 DOS Compile Files Screen	48
Figure 4-28 BoCSharpCompiler Script Finish Screen	49
Figure 4-29 Migration Completed Successfully Screen	49
Figure 4-30 Results of the Migration Screen	50
Figure 4-31 Data Migration Complete Screen	51
Figure 4-32 Readme.txt - Notepad Screen	52

1.0 SCOPE

1.1 IDENTIFICATION

This guide describes Data Migration and the procedures for migrating data from a Block 1 (B1) system to a Block 2 (B2) system.

1.2 DATA MIGRATION OVERVIEW

Data Migration represents the intentional upgrade of a B1 system to a B2 system. Several steps must be followed for the successful migration of data. It is important to note the B2 installation is not intended to automatically upgrade a B1 system. In addition, it is recommended the higher level echelons of care installations be upgraded to the B2 environment prior to the lower levels of care.

Data Migration is the movement of data from one system to another. The movement of data can be from one machine to another; or it is possible to export the data from a B1 system to an external media storage device or network drive, then import the data to the same machine once it has been updated to the B2 system.

The Data Migration Guide must be read and understood prior to beginning the installation.

1.3 DOCUMENT OVERVIEW

This document provides an understanding of the migration of data from a B1 system to a B2 system. It is organized into four sections and two appendixes.

- Section 1 Scope
- Section 2 Pre-Migration Setup
- Section 3 Data Migration Export
- Section 4 Data Migration Import
- Appendix A Acronyms and Abbreviations Table
- Appendix B Pre-Migration Checklist

2.0 PRE-MIGRATION SETUP

Data Migration represents the movement of data located on a current B1 system to an external media source. Once collected, this data can then be migrated to a B2 system. Prior to the migration occurring, the System Administrator (SA) must ensure the following pre-migration items have been followed. A Pre-Migration Checklist has been included in Appendix B.

2.1 DATA TO BE MIGRATED

The Theater Medical Information Program-Joint (TMIP-J) Data Migration Utility will only collect and transfer data found in the My Documents folder of the B1 system and any data contained in the

AHLTA-Theater (AHLTA-T) application database. AHLTA-T application configuration data residing outside the database, Windows security logs, passwords, and accounts external to the database will not be migrated. Only files with .doc, .ppt, .xls, .vsd, .mpp, .dot, and .pdf extensions will be migrated.

The SA should move any documents or files to the My Documents folder or save to an external media source so the data can be migrated. Older documents which are not intended to be migrated should be removed from the My Documents folder. Profile specific information of the owner's My Document folder is lost. All files from all users' My Documents folders are saved to the designated My Document folder.

2.2 B1 DATA TRANSFER VIA THE STEM

On the B1 system, it is important that any data waiting to transmit due to the site being configured as stand-alone or because external communications are down be moved along or archived prior to migrating to B2. To remove Internet Information Server (IIS) unsent Bad Mail messages, the SA will have the Bad Mail Utility run in the background under the account having system administrator's privileges. TMIP-J sent messages will stay in the database. The user must perform a data transfer via the Storage Transfer Encapsulation Method (STEM).

2.2.1 B1 STEM Procedure

When the TMIP-J site is configured as a stand-alone or the external communications network is down, TMIP-J must ensure its clinical data is transferred to the Force Health Protection Interim Theater Database (FHP ITDB). This data transfer is accomplished via STEM.

IMPORTANT WARNING – DO NOT run STEM data transfer procedures from a TMIP-J client configuration. STEM procedures are only applicable to TMIP-J server configurations.

STEM volumes are produced on removable media which replaces e-mail as the means for transfer of encrypted messages to and from other TMIP-J sites. The STEM volumes are physically transported to another TMIP-J site that has access to theater communications. Processing of the clinical messages (e.g., compression, encryption, decryption, decompression, etc.) remains the same regardless if the TMIP-J site is configured as a

stand-alone (no theater communications connectivity) or non-stand-alone (theater communications is available) site.

STEM has two operations – Import and Export. The Import operation transfers message files from the removable media into the TMIP-J workstation's database. The Export operation transfers message files from the TMIP-J workstation's database to the removable media. The **STEM Application** screen is shown in Figure 2-1.

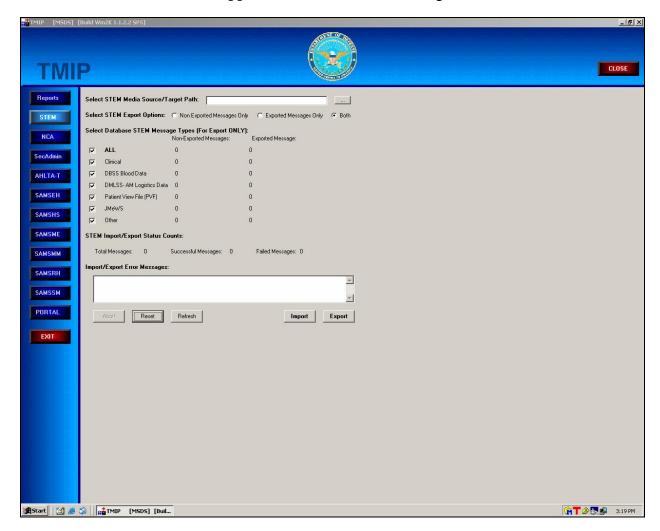


Figure 2-1 STEM Application Screen

2.2.2 Export STEM Volume

This section describes the steps for exporting the message files from the TMIP-J workstation's database to removable media. Refer to Figure 2-1.

- 1. Log into TMIP with a user account having the ROLE_TMIP_STM role.
- 2. Click on the **STEM** button. The **STEM Application** screen is displayed.
- 3. Click on the ... button next to the **Path** field to select the path to the removable media.

- 4. Three options may be selected for exporting to the STEM. Check all of the boxes that apply.
 - When the **Delete on Import/Export** is checked, the messages are deleted from the workstation's database when the export has been successfully completed.
 - When the Re-Export/Import is checked, messages that have previously been exported from the workstation's database but not deleted can be exported again.
 - When the **Abort on Error** is checked, the STEM process halts when an error occurs. The error message displays in the error message box. If this checkbox is not checked, the error displays in the error message box, but processing will continue.
- 5. An informational message is displayed with the following details:
 - The number of messages available for export
 - The number of messages that have been exported
 - The total number of messages
 - The number of messages exported successfully
 - The number of messages that failed export
- 6. Click on the **EXPORT** button to start exporting messages. The following options are available:
 - DELETE Clicking on this button deletes messages that have been exported from the workstation's database. The option may be used when the Delete on Import/Export box is not checked.
 - REPORTS Clicking on this button displays an informational screen reporting the number of previously exported messages and the number of messages waiting to be exported.
 - ABORT Clicking on this button halts the import/export process. The messages already processed will not be backed out.
 - RESET Clicking on this button clears the counters, updates non-exported and exported messages numbers, and clears the error message box.

2.2.3 Import STEM Volume

This section describes the steps for transferring message files from removable media into the TMIP-J workstation's database. Refer to Figure 2-1.

- 1. Log into TMIP-J with a user account having the ROLE_TMIP-J_STM role.
- 2. Click on the **STEM** button. The **STEM Application** screen is displayed.
- 3. Click on the ... button next to the **Path** field to select the path to the removable media.
- 4. Three options may be selected. Check all of the boxes that apply:

- When the **Delete on Import/Export** is checked, the messages are deleted from the workstation's database when the import has been successfully completed.
- When the **Re-Export/Import** is checked, messages that have previously been imported from the workstation's database but not deleted can be imported again.
- When the **Abort on Error** is checked, the STEM process halts when an error occurs. The error message displays in the error message box. If this box is not checked, the error displays in the error message box, but processing will continue.
- 5. An informational message is displayed with the following details:
 - The number of messages available for import
 - The number of messages that have been imported
 - The total number of messages
 - The number of messages imported successfully
 - The number of messages that failed import
- 6. Click on the **IMPORT** button to start importing messages. The following options are available:
 - DELETE Clicking on this button deletes messages that have been imported from the removable media. The option may be used when the Delete on Import/Export box is not checked.
 - REPORTS Clicking on this button displays an informational window reporting the number of previously imported messages and the number of messages waiting to be imported.
 - ABORT Clicking on this button halts the import/export process. The messages that have already processed will not be backed out.
 - RESET Clicking on this button clears the counters, updates non-imported and imported messages numbers, and clears the error message box.

2.3 PREPARING A B1 SYSTEM

There are several steps to complete prior to exporting the data from a B1 system. The B1 system must be either a stand-alone or a server. All clients should be closed before Data Migration starts. The migration processes shuts down the database, effectively disconnecting all clients. Data Migration Utility should be run during a scheduled maintenance downtime and NOT during normal operational hours.

2.3.1 Back Up a B1 System

This is simply a precautionary step to protect the data should the system need to be restored. Refer to the B1 System Administration Manual (SAM) for details on how to back up the B1 System.

In addition, if there are any other files or data that are important, it is recommended the general file backup procedures be followed as appropriate.

2.3.2 B1 TMIP-J and Oracle Version

For the Data Migration Utility to complete the migration process, the TMIP-J version being utilized on the B1 system must be Version 1.1.2.2 Service Pack 5 and the Oracle Database must be Oracle 8i Version 8.1.7.4.17. If these versions are earlier versions, the Data Migration will abort and the database will be reset to the original data. As new Service Packs become available, the users will be notified of the compatibility of the Service Pack to the Data Migration Utility. The SA must upgrade the B1 system to the correct versions of TMIP-J and Oracle prior to running the Data Migration Utility.

2.3.3 Pre-Export Migration Setup

Prior to the start of the Data Migration Utility, the SA should create one working folder and one container folder for the export portion of the migration. These folders should be located on an external media source, or the Working folder can be created on a local drive, and the TMIP-J B2 Data Container folder can be created on a network drive. The Working folder pulls the information from the applications and then transfers the data to the TMIP-J B2 Data Container folder where the information is encrypted and then utilized during the import portion of the Data Migration process.

2.4 Preparing a B2 System

2.4.1 Set Up a B2 Machine

Prior to importing B1 data, the B2 machine must be:

- Formatted or imaged
- Hard drives partitioned
- TMIP-J B2 Installation Wrapper installed

The minimum hardware requirements for the common scenario in which a laptop must be able to operate as either a client or as a stand-alone machine are as follows:

- Processor: Pentium III, 800 Megahertz (MHz)
- RAM: 1 Gigabyte (GB)
- Hard Disk: 40 GB of free space

This can be done prior to starting the Data Migration Utility if a different machine is being used as the B2 machine or just prior to importing the data to the B2 machine.

Note: Data Migration Utility will run on any machine with Windows 2003 or XP.

Note: Data Migration Utility will update from a Block 1 Service Pack 5 (B1 SP5) machine.

Refer to the TMIP-J SAM Appendix B: System Installation Guide (SIG) for the installation and setup of a new machine.

2.4.2 Pre-Import Migration Setup

Prior to the start of the Data Migration Utility, the SA should create two folders for the import portion of the migration. These folders should be located on an external media source or a local drive. The folders should be named "Working" and "Restore." The Working folder will pull the information from TMIP-J B2 Data Container folder created in the export portion of the Data Migration, decrypt the information, and then import into the B2 machine. The My Documents folder migrated from the B1 machine will be migrated to the Restore folder. After the migration process is complete, the Working folder will automatically be deleted. A pre-migration checklist is provided in Appendix B.

Note: To retain any patient information or records, AHLTA-T and Shipboard Non-Tactical Automated Data Processing Program Automated Medical System (SAMS) should not be used during the Data Migration process. Since the Data Migration Utility will overwrite any information in either application, users should not input any data in either AHLTA-T or SAMS prior to the import of data.

2.5 POST-DEPLOYMENT SOFTWARE SUPPORT

The Military Health System (MHS) provides Service Help Desk operational sustainment support. In general, any software issues should be reported first to the Service Help Desk. The Service Help Desk will then escalate issues up through the MHS Help Desk tier levels as deemed necessary.

TMIP-J supports the Service Help Desk including operational sustainment support and mid-cycle, low risk and maintenance releases for technology updates and problem fixes. The site SA will provide the Service a Help Desk telephone number once the software is deployed. The MHS Service Help Desk is available via the following:

- http://www.mhs-helpdesk.com
- help@mhs-helpdesk.com
- 800.600.9332 (Continental United States [CONUS])
- **210.767.5250 (Direct)**
- 866.637.8725 (Outside Continental United States [OCONUS])

3.0 DATA MIGRATION – EXPORT

The first step in Data Migration is the exporting of data from the B1 system. The primary migration strategy is to export all relevant patient data from the B1 AHLTA-T system. The following outlines the procedures necessary to export data from the B1 AHLTA-T system.

The Data Migration Utility expected runtime for the exporting of a small database (0-10,000 encounters) is approximately 5 minutes, a medium sized database (11,000-25,000 encounters) is approximately 10 minutes, and a large database (26,000-50,000 encounters) is approximately 20 minutes.

Note: The TMIP-J version needs to be at least Version 1.1.2.2 Service Pack 5 (SP5) for the migration process to work.

Note: During the migration process, all users should log off either AHLTA-T or SAMS applications and close the applications. The Data Migration Utility should be run during scheduled maintenance downtime.

Note: The Data Migration process could fail if the user manually edits any entry in the tnsnames.ora file that has a SERVICE_NAME = ldb. Data Migration may not be able to properly determine if the database is configured as a remote client. The user should only use the provided tools to modify system configuration. On a B1 system, use the TMIP Configuration tool and on a B2 system, use the Oracle TnsNames Manager tool. See the TMIP-J SAM for additional information.

- 1. Use Windows Explorer from Windows Desktop to navigate to the CD/DVD drive.
- 2. Double click on **setup.exe** file of the Data Migration Utility (**Figure 3-1**). Data Migration Utility is found on B2FY07SP1-DVD04 DVD.

Note: Data Migration Utility **CANNOT** be executed from a network location.

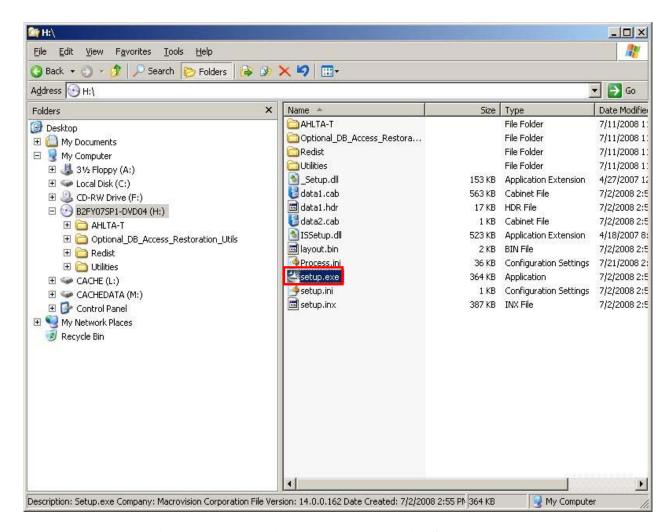


Figure 3-1 Data Migration setup.exe File Screen

- The Data Migration Utility will detect the current TMIP-J version.
 - If the TMIP-J version is 1.1.2.2, the Data Migration Utility will continue.
 - If the TMIP-J version is not 1.1.2.2, the Data Migration Utility will ask the user to upgrade to TMIP-J version 1.1.2.2.
 - If the Data Migration Utility can not detect TMIP-J version, the process will abort.

Note: .NET Framework only needs to be installed the first time the Data Migration Utility is run if the system does not already have .NET Framework installed.

Note: Following screens will appear during .Net Framework installation.

- 3. Microsoft .NET Framework 1.1 Setup License Agreement screen appears (Figure 3-2).
- 4. Select **I agree** radio button.
- 5. Click Install.

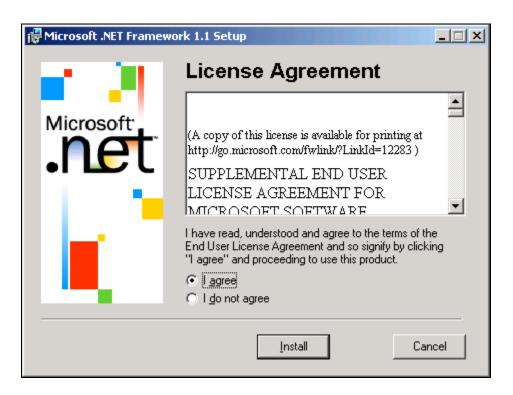


Figure 3-2 .NET Framework 1.1 Setup License Agreement Screen

- 6. **Installation of Microsoft .NET Framework 1.1 is complete** screen appears (**Figure 3-3**).
- 7. Click OK.



Figure 3-3 Installation of Microsoft .NET Framework 1.1 Is Complete Screen

Step

- 8. Welcome to the Data Migration Wizard screen appears (Figure 3-4).
- 9. Select Next.

WARNING: Attempting to run Data Migration on a system more than once will cause a dialog screen to appear stating, Continuing with this application may corrupt databases and/or other files. Do you want to continue?



Figure 3-4 Welcome to the Data Migration Wizard Screen

Step 10. Data Migration Guide dialog box appears (Figure 3-5). 11. Select OK.



Figure 3-5 Data Migration Guide Dialog Box

Step 12. Data Migration remote users warning screen appears (Figure 3-6). 13. Select OK.



Figure 3-6 Data Migration Remote Users Warning Screen

Step

14. Choose Working And TMIP B2 Data Container Folders screen appears (Figure 3-7).

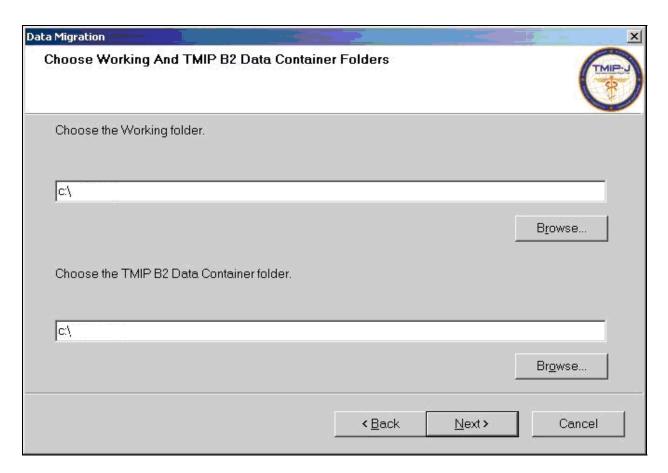


Figure 3-7 Choose Working And TMIP B2 Data Container Folders Screen

- 15. Enter Working in the Choose the Working folder field.
- 16. Enter Container in the Choose the TMIP B2 Data Container folder field.
- 17. Click Next.

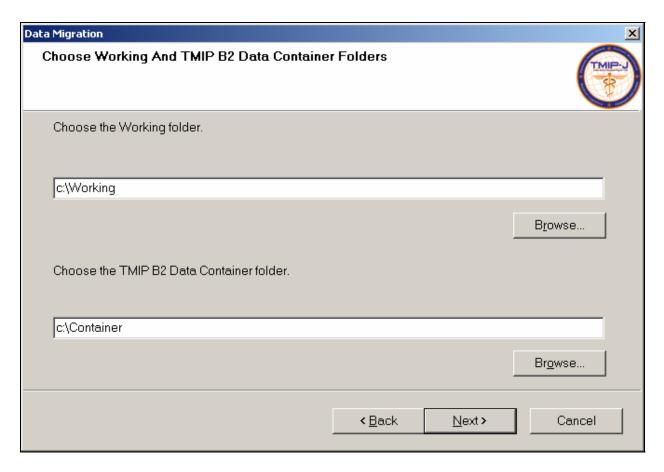


Figure 3-8 Choose Working And TMIP B2 Data Container Folders - Completed Screen

- 18. **Data Migration** drive selection screen appears (**Figure 3-9**).
- 19. If **Yes**, a Network Drive can be chosen.
- 20. If No, Data Migration User input screen appears (Figure 3-10).



Figure 3-9 Data Migration Drive Selection Screen

21. Enter Oracle System UserID - system. 22. Enter Oracle System Password. 23. Enter SID - ldb. 24. Click OK.

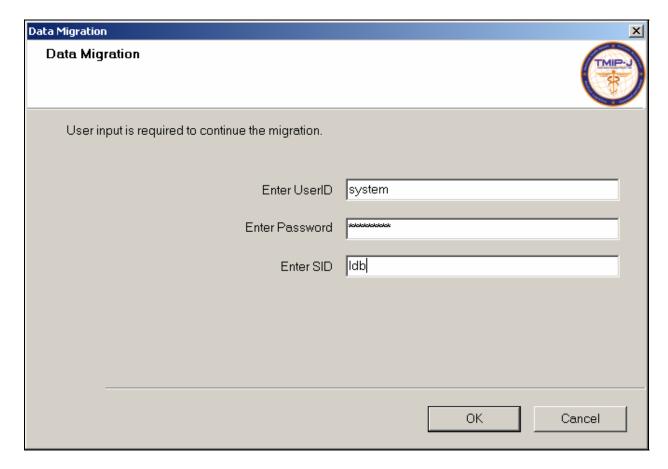


Figure 3-10 Data Migration User Input Screen

- 25. Select the Product to migrate from: screen appears (Figure 3-11).
- 26. Select TMIP B1 SP5.
- 27. Click Next.

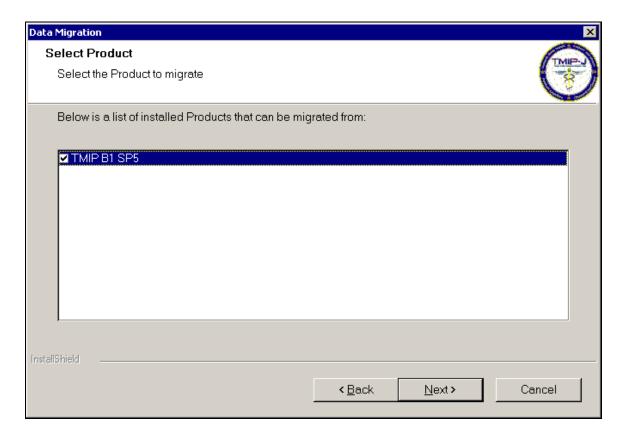


Figure 3-11 Select the Product to Migrate Screen

28. List of Components to migrate screen appears (Figure 3-12). 29. Click Next.

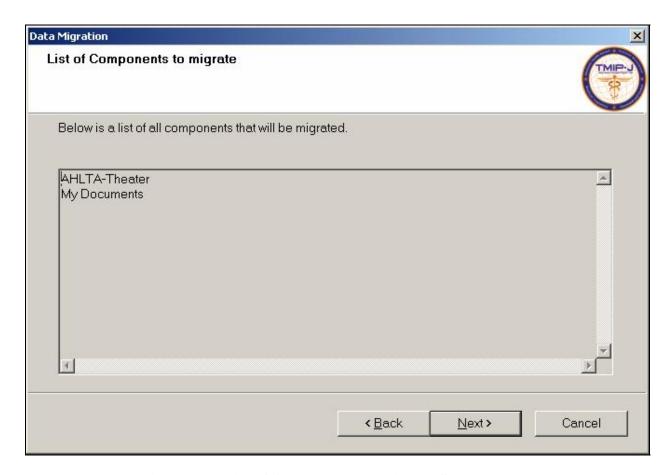


Figure 3-12 List of Components to Migrate Screen

Step 30. Data Migration Setup Status screen appears (Figure 3-13).

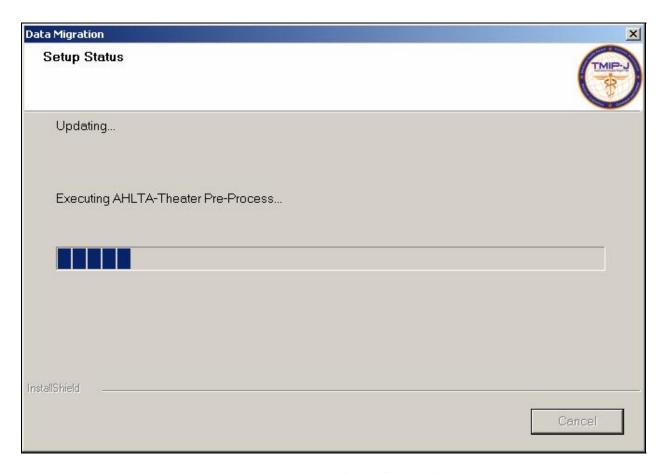


Figure 3-13 Data Migration Setup Status Screen

- 31. The Data Migration process will require the approximate amount of time to complete screen appears (Figure 3-14).
- 32. Click **OK**.

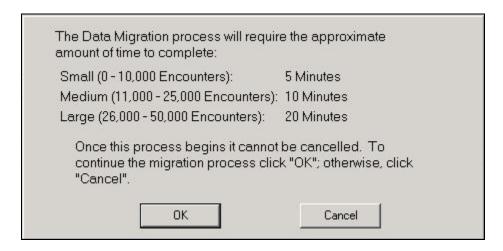


Figure 3-14 Data Migration Estimated Time to Complete Screen

Step

33. **CHCSII-Theater Data Migration Manager** screen appears (**Figure 3-15**) showing number of patches completed and the estimated time remaining for export.

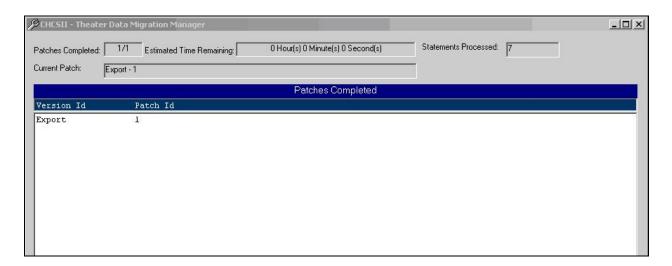


Figure 3-15 CHCSII-Theater Data Migration Manager Screen

34. Oracle DOS Export Data screen for exporting data launches (Figure 3-16).

```
🚾 D:\oracle\ora81\bin\exp.exe
                                                                                                                                       definitions for user GEMS
definitions for user GEMSA
definitions for user C2T_X
definitions for user SCADB
definitions for user SECADMIN
    exporting object
                                 type
   exporting object exporting object
                                 type
   exporting object type exporting object type
exporting object type definitions
exporting object type definitions
About to export GEMS's objects ...
exporting database links
exporting sequence numbers
exporting cluster definitions
about to export GEMS's tables via
exporting table
                                                   les via Conventional Path ...
ABNORMAL_CODES
ACCESS_AUDIT
ACUTE_DIAGNOSES_CODES
                                                                                                         84 rows
                                                                                                                      exported
       exporting exporting
                         table
                                                                                                              rows
                                                                                                                       exported
                         table
                                                                                                              rows
                                                                                                                       exported
       exporting
                                                                                                                       exported
                                                                                                              rows
       exporting
                                                                                                              rows
                                                                                                                       exported
       exporting
                                                ALERT_ADE_DOCEUT
                                                                                                              rows
                                                                                                                       exported
       exporting
                                                                                                              rows
                                                                                                                       exported
       exporting
                                                                                                               rows
       exporting
                                                                                                               rows
       exporting
                                                                                                               rows
                                                                                                                       exported
       exporting
                                                                                                              rows
                                                                                                                       exported
       exporting exporting
                                                                ALLERGY COMMENT
                         table
                                                                                                              rows
                                                                                                                       exported
                                                     ALLERGYCOMMENT_AUDIT
                         table
                                                                                                              rows
                                                                                                                       exported
                                                              ALLERGY REACTION
       exporting
                                                                                                                      exported
                         table
                                                                                                              rows
                                                    ALLERGYREACTION_AUDIT
       exporting
```

Figure 3-16 Oracle DOS Export Data Screen

35. CHCS II-Theater process completed successfully screen appears (Figure 3-17). 36. Click OK.



Figure 3-17 CHCS II-Theater Process Completed Successfully Screen

37. Various sequence of **Data Migration Setup Status** screens appear (**Figure 3-18**) and (**Figure 3-19**).

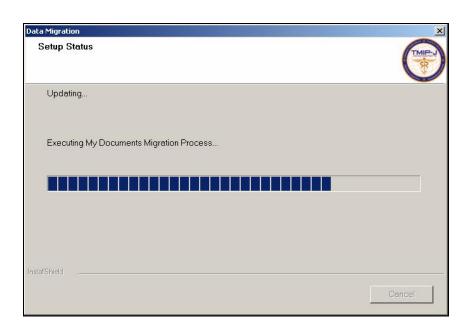


Figure 3-18 Data Migration Setup Status Screen

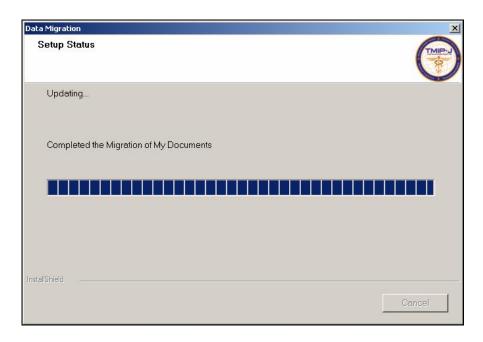


Figure 3-19 Data Migration Setup Status Screen - Continued

38. Task Status screen appear displaying status of Oracle backup (Figure 3-20).

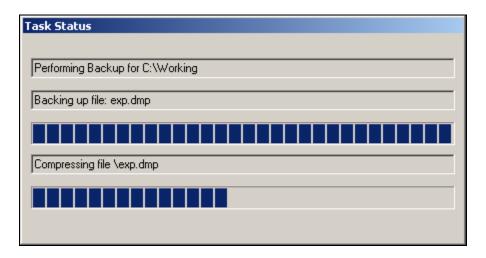


Figure 3-20 Task Status Screen

- 39. The Migration completed successfully screen appears (Figure 3-21).
- 40. Click **OK**.

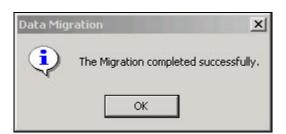


Figure 3-21 Migration Completed Successfully Screen

- 41. **Results of the Migration** screen opens indicating Data Migration successes or failures (**Figure 3-22**).
- 42. Click Next.

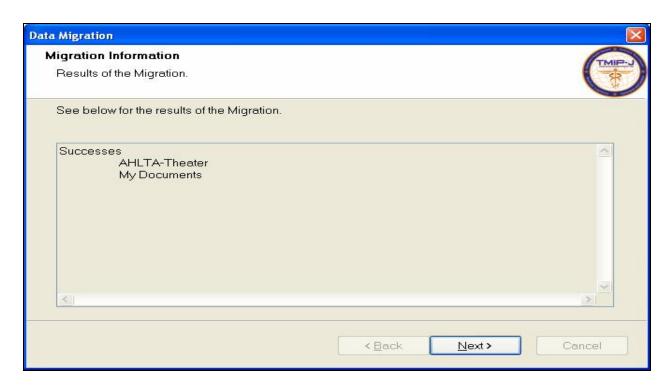


Figure 3-22 Results of Migration Screen

- 43. Data Migration Complete screen appears (Figure 3-23).
- 44. Click Finish.

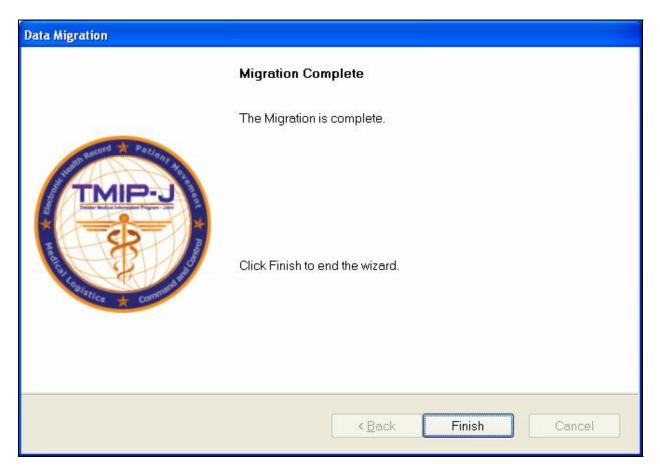


Figure 3-23 Migration Complete Screen

Step

45. Close Readme.txt - Notepad screen (Figure 3-24).

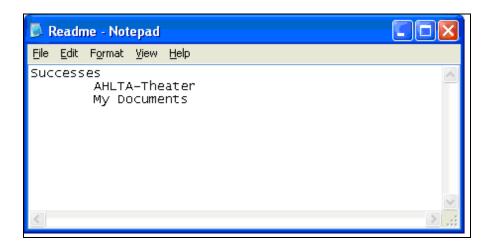


Figure 3-24 Readme.txt Notepad Screen

4.0 DATA MIGRATION – IMPORT

After a successful export of data from the B1 system, the SA must then import the data into the new B2 machine. The encrypted exported data is decrypted and migrated to the new system. Any information saved and migrated in the My Documents folder will be saved to the Restore folder created prior to running the Data Migration Utility.

Data Migration Utility expected runtime for importing of a small database (0-10,000 encounters) is approximately 1 hour and 35 minutes, a medium sized database (11,000-25,000 encounters) is approximately 2 hours and 15 minutes, and a large database (26,000-50,000 encounters) is approximately 3 hours and 15 minutes. If the Data Migration Utility has been successfully run on the B2 machine, the SA receives a warning window stating, "Migration has been completed on this machine. Do you want to continue?" If the SA chooses to continue, the data on the machine is overwritten with the new data.

Note: The TMIP-J version needs to be at least Version 1.1.2.2 Block 2 Release 1 Service Pack 2 (B2 R1 SP2) for the migration process to work.

Note: AHLTA-T and SAMS should not be used prior to the import portion of the Data Migration Utility because any data entered into either application will be overwritten during the import process on the B2 System.

Note: Data Migration process could fail if the user manually edits any entry in the tnsnames.ora file that has a SERVICE_NAME = ldb. Data Migration may not be able to properly determine if the database is configured as a remote client. The user should only use the provided tools to modify system configuration. On a B1 system, use the TMIP-J Configuration tool and on a B2 system, use the Oracle TnsNames Manager tool. See the TMIP-J SAM for additional information.

- 1. Use Windows Explorer from Windows Desktop to navigate to the CD/DVD drive.
- 2. Double click on **setup.exe** file of Data Migration Utility (**Figure 4-1**). Data Migration Utility is found on B2FY07SP1-DVD04 DVD.

Note: Data Migration Utility **CANNOT** be executed from a network location.

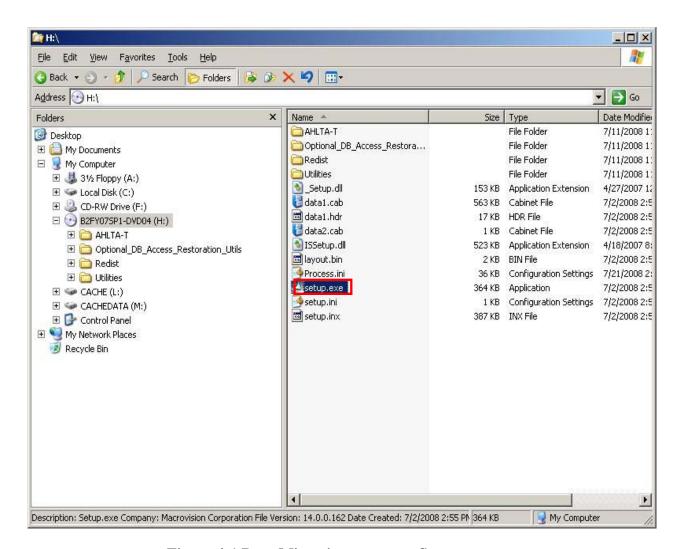


Figure 4-1 Data Migration setup.exe Screen

- 3. Welcome to the Data Migration Wizard screen appears (Figure 4-2).
- 4. Select Next.



Figure 4-2 Welcome to the Data Migration Wizard Screen

Step

Note: If the following **Error** message screen appears (**Figure 4-3**), migration process cannot continue because (TMIP B1 SP5 or TMIP B2 R1 SP2 or TMIP B2 FY07 SP1) product is not installed.

By selecting **OK**, the install process will terminate.

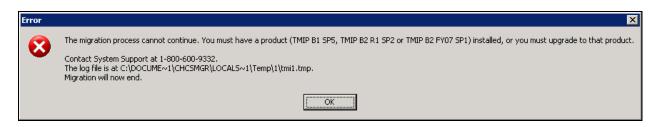


Figure 4-3 Data Migration Error Screen

- 5. Data Migration Guide dialog screen appears (Figure 4-4).
- 6. Select OK.



Figure 4-4 Data Migration Guide Dialog Screen

- 7. **Data Migration remote users warning** screen appears (**Figure 4-5**).
- 8. Select OK.

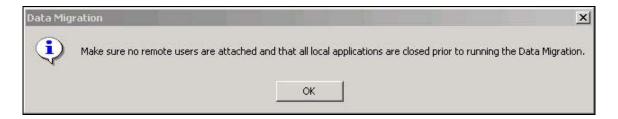


Figure 4-5 Data Migration Remote Users Warning Screen

- 9. Choose Working Folder, MyDocs Restore Folder And TMIP B2 screen appears (Figure 4-6).
- 10. Select the browse button [...] to locate **TMIP B2 Data Container file.**

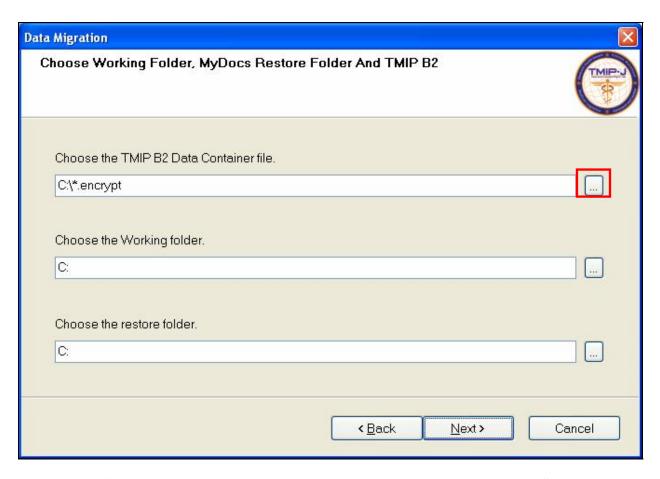


Figure 4-6 Choose Working Folder, MyDocs Restore Folder And TMIP B2 Screen

SAIC D/SIDDOMS Doc. D3-TMIP-5019A 20 Sep 2008

- 11. Select **TMIPB2DataContainer** folder from the **Select File** screen. This will prefill the **Choose the TMIP B2 Data Container file** field (**Figure 4-7**).
- 12. Enter WorkingFolder in the Choose the Working folder field.
- 13. Enter RestoreFolder in the Choose the restore folder field.
- 14. Click Next.

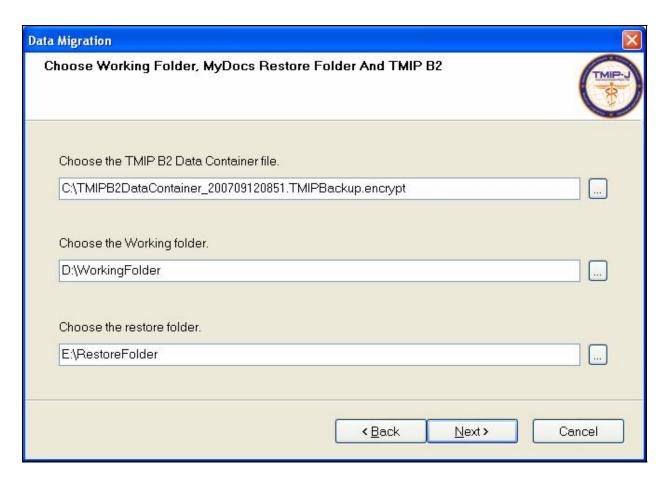


Figure 4-7 Select TMIPB2DataContainer Folder from Select File – Completed Screen

Note: Depending on which version of TMIP-J software running, either Figure 4-8 or Figure 4-9 will display. Follow the step listed under the correct version displayed figure.

- 15. Select the Product to migrate to: screen appears (Figure 4-8).
- 16. Select TMIP B2 R1 SP2.
- 17. Click Next.

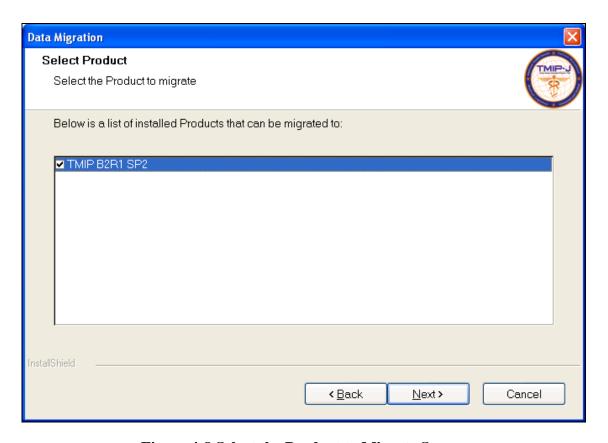


Figure 4-8 Select the Product to Migrate Screen

18. List of Components to migrate screen appears (Figure 4-11). 19. Click Next.

- 20. **Select the Product to migrate to**: screen appears (Figure 4-9).
- 21. Select TMIP B2 FY07 SP1.
- 22. Click Next.

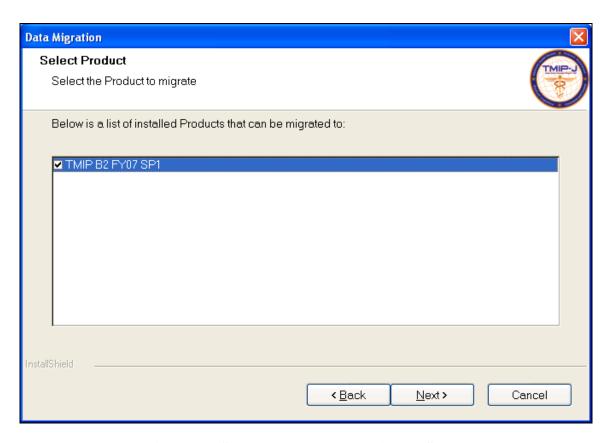


Figure 4-9 Select the Product to Migrate Screen

- 23. List of Components to migrate screen appears (Error! Reference source not found.).
- 24. Click Next.

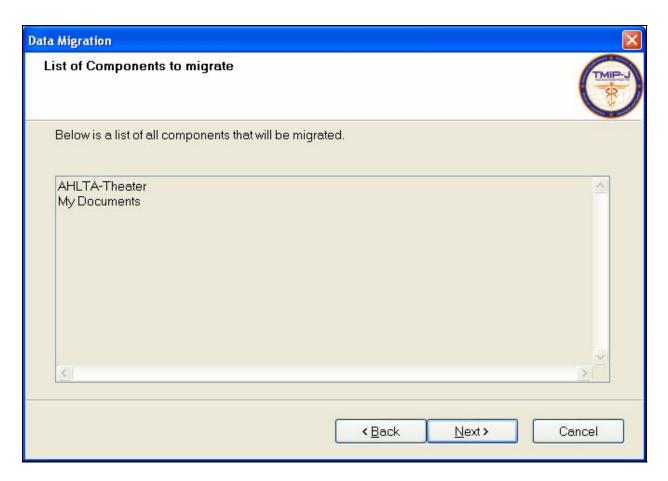


Figure 4-10 List of Components to Migrate Screen

SAIC D/SIDDOMS Doc. D3-TMIP-5019A 20 Sep 2008

- 25. Data Migration User input screen appears (Figure 4-11).
- 26. Enter Oracle System UserID system.
- 27. Enter Oracle System Password.
- 28. Enter SID ldb.
- 29. Click **OK**.

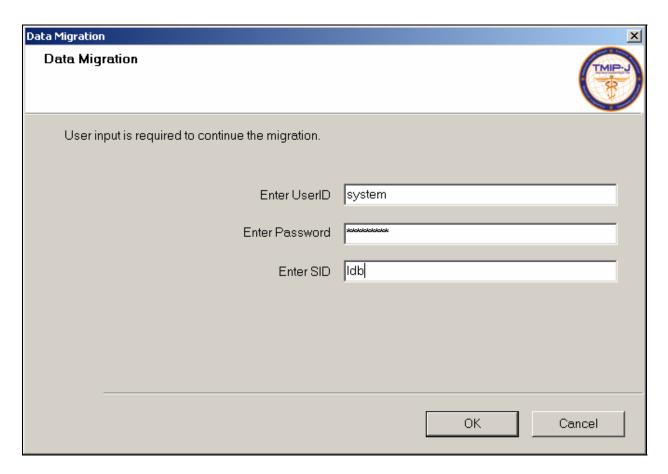


Figure 4-11 Data Migration User Input Screen

- 30. Task Status screens showing status of Oracle backup appears (Figure 4-12) and (30
- 31. **Figure 4-13**). No user interaction for these screens is required.


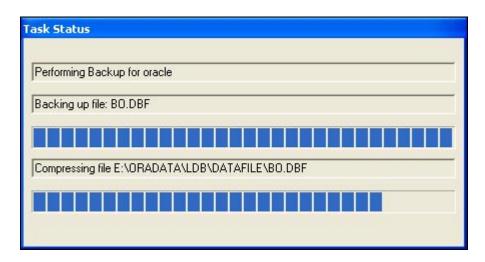


Figure 4-12 Task Status Screen

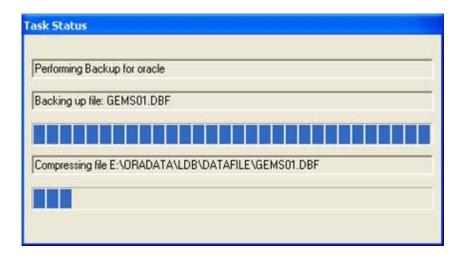

30

Figure 4-13 Task Status Screen - Continued

32. **Setup Status** screen displays with status bar advancing as the data is imported (**Figure 4-14**).

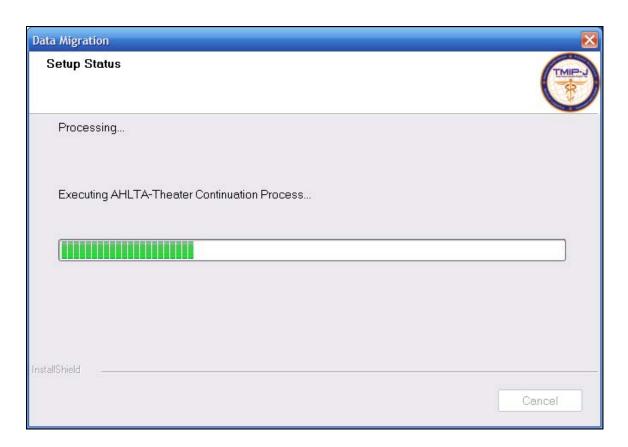


Figure 4-14 Setup Status Screen

SAIC D/SIDDOMS Doc. D3-TMIP-5019A 20 Sep 2008

Step

- 33. Existing data in the application database warning dialog box appears (Figure 4-15).
- 34. Click **OK**.

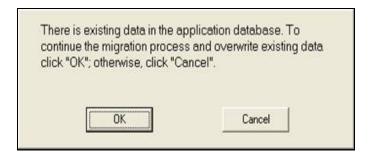


Figure 4-15 There Is Existing Data in the Application Database Warning Dialog Box

- 35. Data Migration estimated time to complete screen launches (**Figure 4-16**).
- 36. Click **OK**.

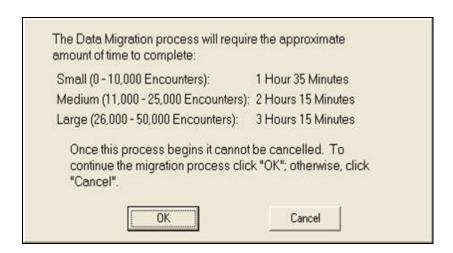


Figure 4-16 Data Migration Estimated Time to Complete Screen

37. **AHLTA-Theater Data Migration Manager** screen (**Figure 4-17**) opens to show number of patches completed and the estimated time remaining for import.

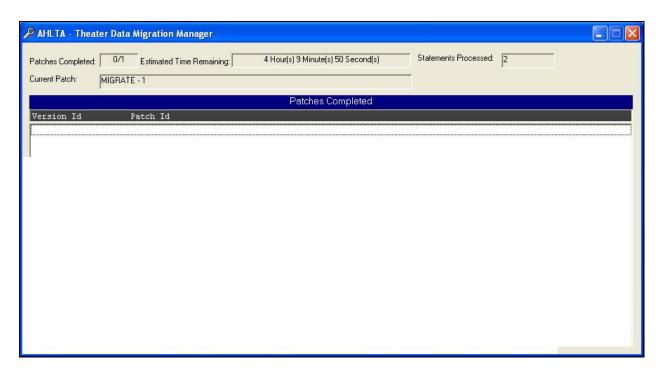


Figure 4-17 AHLTA-Theater Data Migration Manager Screen

38. DOS window for Oracle Product Import screen appears (**Figure 4-18**).



Figure 4-18 DOS Screen for Oracle Product Importing

SAIC D/SIDDOMS Doc. D3-TMIP-5019A 20 Sep 2008

Step

39. **AHLTA-Theater Data Migration Manager** screen (**Figure 4-19**) opens to show number of patches completed and the estimated time remaining for import.


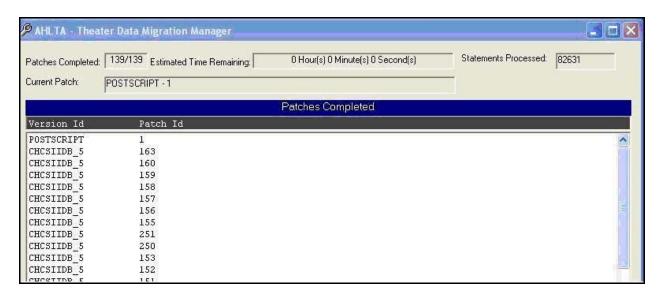


Figure 4-19 AHLTA-Theater Data Migration Manager Screen

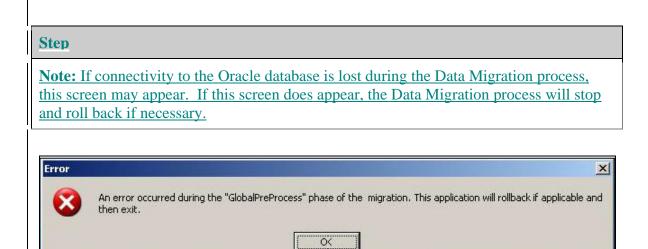

- 40. AHLTA-Theater process completed successfully screen appears (Figure 4-20).
- 41. Click OK.



Figure 4-20 AHLTA-Theater Process Completed Successfully Screen

42. The following screens are a sequence of setup status screens that require no user interaction

(Figure 4-21 through Figure 4-24).

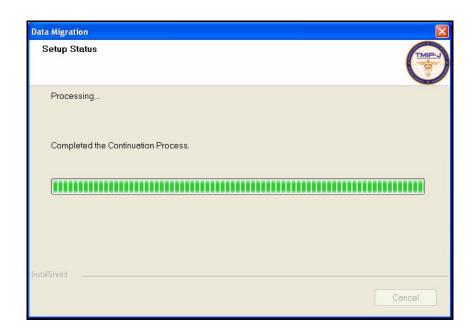


Figure 4-21 Setup Status Continuation Process Screen

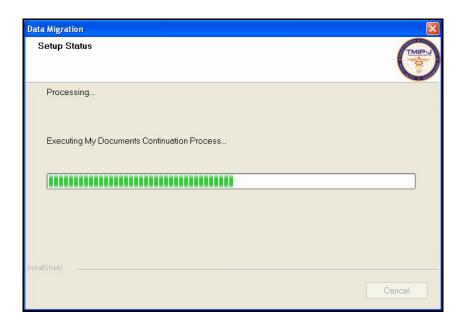


Figure 4-22 Setup Status My Documents Continuation Process Screen - Continued

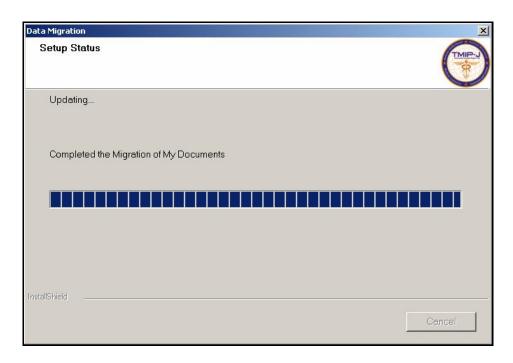


Figure 4-23 Completed Migration of My Documents Screen - Continued

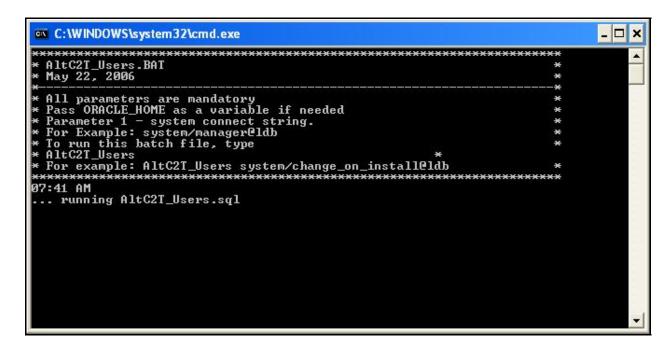


Figure 4-24 DOS Compile Screen

Step
43. BoCSharpCompiler Please Logon screen appears (Figure 4-25).

44. Enter Sys Admin User ID - sys. 45. Enter Password. 46. Enter SID - LDB. 47. Click OK.

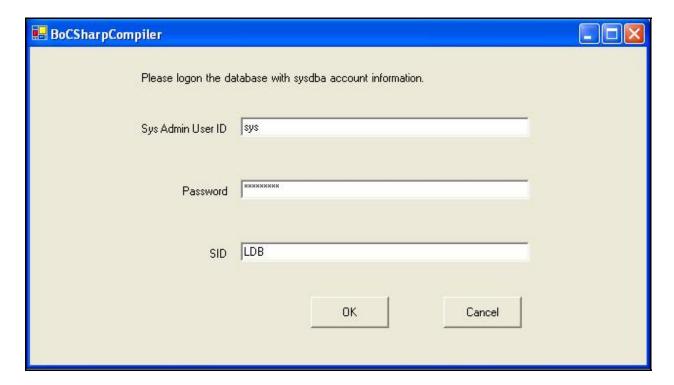


Figure 4-25 BoCSharpCompiler Logon Screen

48. **BoCSharpCompiler logon** to database screens appear (**Figure 4-26**) and (**Figure 4-27**).

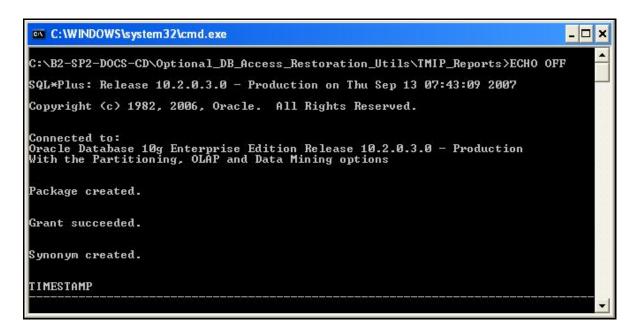


Figure 4-26 BoCSharpCompiler Package Screen

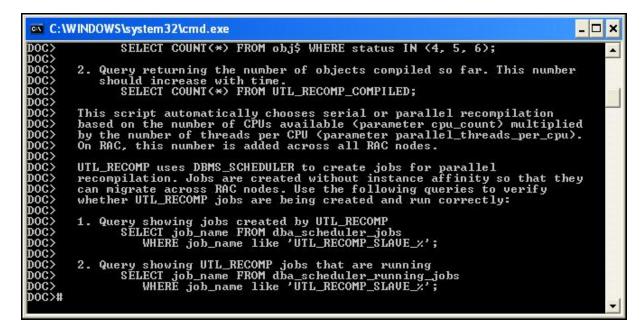


Figure 4-27 DOS Compile Files Screen

- 49. BoCSharpCompiler Script Finish screen appears (Figure 4-28).
- 50. Click **OK**.

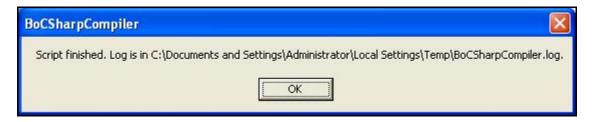


Figure 4-28 BoCSharpCompiler Script Finish Screen

- 51. The Migration completed successfully screen appears (Figure 4-29).
- 52. Click **OK.**



Figure 4-29 Migration Completed Successfully Screen

- 53. **Results of the Migration** screen (**Figure 4-30**) displaying data migration successes or failures appear.
- 54. Click Next.

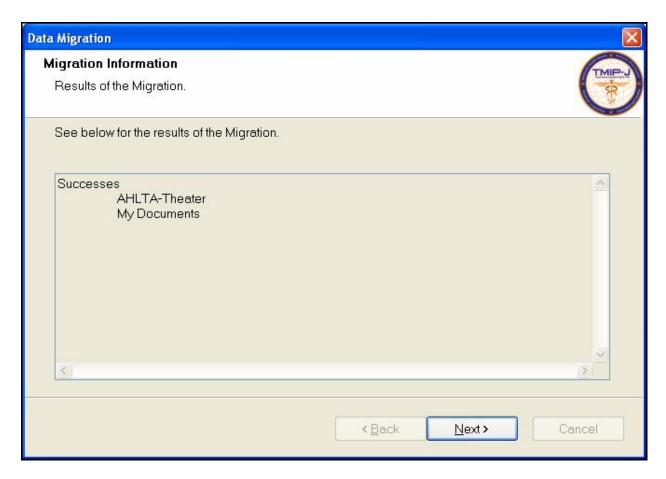


Figure 4-30 Results of the Migration Screen

- 55. Data Migration Complete screen appears (Figure 4-31).
- 56. Click Finish.



Figure 4-31 Data Migration Complete Screen

57. Close **Readme.txt** - **Notepad** screen (**Figure 4-32**).

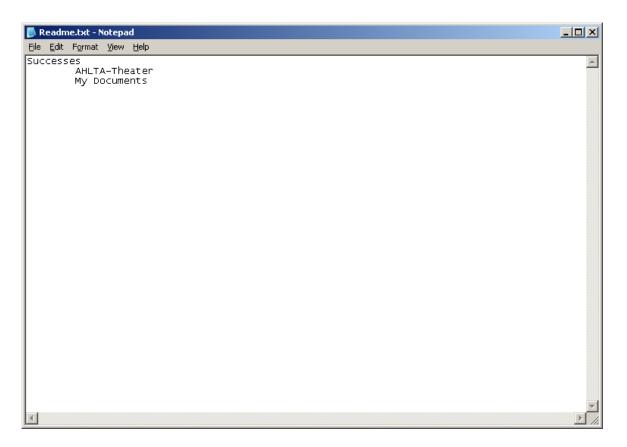


Figure 4-32 Readme.txt - Notepad Screen

Step

- 58. After successful migration, all files from the Data Migration working folder, excluding logs, will be removed as part of post clean up process.
- 59. If migration fails, the database will be restored to the previous backup data.

After a successful Data Migration install, the System Administrator needs to log into Data Migration. During the Data Migration install, the System Administrator's password is automatically set to the default password. It is this password, the default password, that the System Administrator needs to use to log into the system.

APPENDIX A: ACRONYMS/ABBREVIATIONS

Acronym/Abbreviation	Definition		
AHLTA-T	AHLTA-Theater		
B1	Block 1		
B2	Block 2		
B2 FY07 SP1	Block 2 Fiscal Year 07 Service Pack 1		
B2 R1 SP2	Block 2 Release 1 Service Pack 2		
CONUS	Continental United States		
FHP ITDB	Force Health Protection Interim Theater Database		
GB	Gigabyte		
IIS	Internet Information Server		
MHS	Military Health System		
MHz	Megahertz		
OCONUS	Outside Continental United States		
PMO	Program Management Office		
SA	System Administrator		
SAM	System Administration Manual		
SAMS	Shipboard Non-Tactical Automated Data Processing Program Automated Medical System		
SID	Security Identifier		
SIG	System Installation Guide		
SP5	Service Pack 5		
STEM	Storage Transfer Encapsulation Method		
TMIP-J	Theater Medical Information Program-Joint		

APPENDIX B: PRE-MIGRATION CHECKLIST

Item	Task	Reference Section	Complete
1	My Documents folder contains all data to be migrated.	2.1	
2	The IIS Bad Mail folder and COMMS table are empty.	2.2	
3	Prepare the B1 machine is to be upgraded to a B2 machine.	2.3	
4	B1 machine configured as a stand-alone.	2.3	
5	Perform system backup of the B1 system.	2.3.1	
6	B1 machine has TMIP-J version 1.1.2.2 Service Pack 5 or later version.	2.3.2	
7	Oracle 8i version 8.1.7.4.17 or later on the B1 system.	2.3.2	
8	Create the Working folder for exporting data.	2.3.3	
9	Create the TMIP-J B2 Data Container folder for exporting data.	2.3.3	
10	External media source identified for the Data Migration.	2.3.3/2.4.2	
11	Prepare B2 Machine.	2.4.1	
12	Create the Working folder for importing data.	2.4.2	
13	Create the Restore folder for importing data.	2.4.2	
14	Log off all remote users on the AHLTA-T and SAMS applications and close all applications.	3	
15	Begin Data Migration Export of B2 data.	3	